

## Axiom™ Patient Care Headwall

The Axiom™ Patient Care Headwall is manufactured in a variety of configurations and sizes to meet the specific needs of the various areas in the hospital.

### **The Half-wall**

The Half-wall is extremely useful in a room where a concentration of patient services [electrical, medical gas, nurse call, telephone] is needed and only limited space is available. It can service one or two patients. It consists of a single section and mounts approximately 30 to 45 inches above the finished floor.

### **The Single Section**

The Single section headwall can be used in a general or critical care area. It can be equipped with a wide variety of services, including electrical, medical gas, communications, monitor, etc. It can also service one or two patients. The unit pictured is equipped to service one patient. The single section headwall is the most popular version of the various headwall systems because it offers great flexibility in selection of patient care services to meet a wide range of applications.



### **Double Section Units**

The double section headwall consists of two vertical sections mounted adjacent to each side of the bed. It is often used for one patient when the single section does not offer enough space or accessibility for patient services. The double section headwall permits a greater selection of patient services, and allows splitting these services to serve the patient from either side – allowing for more efficient use by hospital personnel.

### **Triple Section Units**

The triple section headwall consist of three separate sections containing a variety of patient services.


©HSI 2009D-2155

### Features

- Standard vertical equipment tracks are part of each headwall
- Large selection of track accessories
- Flexible headwall designs – easily accepts customer options
- Durable, low maintenance construction
- Accepts all brands of medical gas outlets or nurse call systems

### Options

- Nurse call
- Medical Gas Outlets
- Equipotential Grounding system
- Storage cabinets
- Bed bumpers
- Monitor support
- Isolated Power
- Emergency Electrical Receptacles
- Normal Electrical Receptacles
- And many more

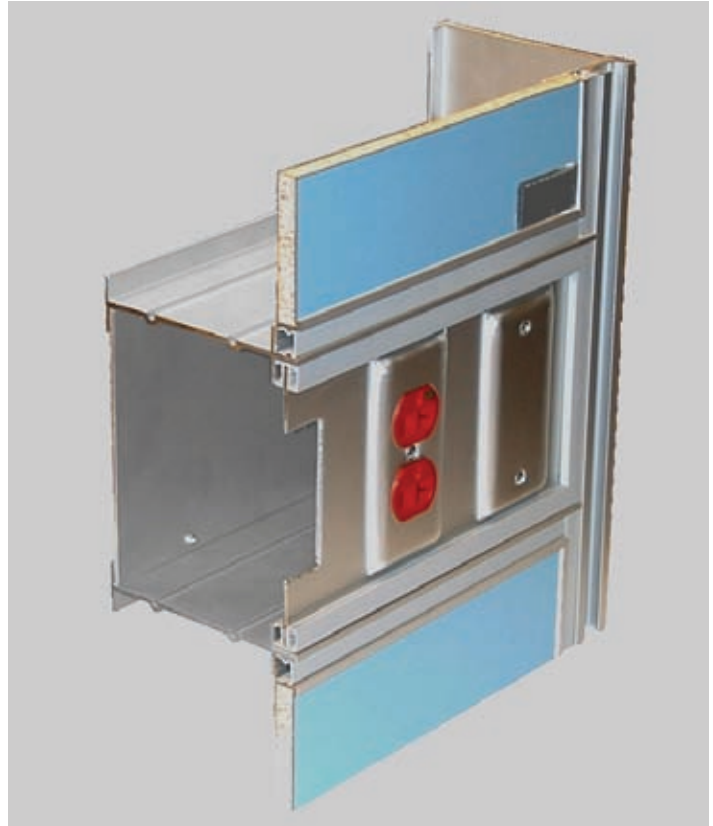

OSHPD  
R-0136